

# USE AND CARE

Congratulations on your purchase of the Obus Ultra Forme™ Seat. The Obus Ultra Forme™ Seat's unique design evenly distributes your body weight, assists circulation through the lower extremities and reduces pressure points by adding cushioning comfort to your tailbone, sitting bones and thighs.

**When used with the Obus Ultra Forme™ Backrest, it will create a complete, comfortable and ergonomic sitting experience by:**

- Encouraging optimum alignment of the pelvis and spine
- Reducing pressure and postural strain on the back

## FEATURES


- Contoured sides support your pelvis and hips in a comfortably balanced position.
- Raised center portion separates your legs slightly and tilts your pelvis backward to eliminate excess lumbar curve and reduce spinal strain.
- Cut-out at the center of the seat's posterior reduces pressure on your tailbone.
- Removable front zippered section adjusts the depth to fit almost any chair, allowing for proper knee clearance and unrestricted circulation.
- Two "hook and loop" fastening tabs allow you to securely attach your Obus Ultra Forme™ Seat to the Obus Ultra Forme™ Backrest or any Obus Forme™ by Conair® Backrest Support.
- Polyurethane foam construction absorbs compression forces and dissipates the vibrations created in moving vehicles, making it the perfect traveling companion.

## CREATING YOUR ERGONOMIC SEATING SYSTEM

Create a complete, comfortable and ergonomic sitting experience by using your Obus Ultra Forme™ Seat together with the Obus Ultra Forme™ Backrest. When used together, they encourage optimum alignment of the pelvis and spine, and reduce pressure and postural strain on the back.

## USING YOUR ERGONOMIC SEATING SYSTEM

1. To attach the Obus Ultra Forme™ Seat to the Obus Ultra Forme™ Backrest, take the two "hook and loop" fastening tabs located at the back of the seat and position and press them up and onto the back of the Obus Ultra Forme™ Backrest for a snug fit (see diagrams).
2. Move your ergonomic seating system into its location and position it with the Obus Ultra Forme™ Backrest resting on top of the Obus Ultra Forme™ Seat.
3. Stretch the elastic straps of the Obus Ultra Forme™ Backrest around the chair or seat. Tighten the straps using the adjustable buckles. Again, ensure the backrest is resting on top of the seat.
4. Adjust the position of the back of your chair by tilting it backward or forward as necessary.
5. Before using your Obus Ultra Forme™ Seat or the Obus Ultra Forme™ Backrest in a vehicle, please consult the owner's manual to ensure that you can properly adjust your vehicle's seat and/or headrest after installation.

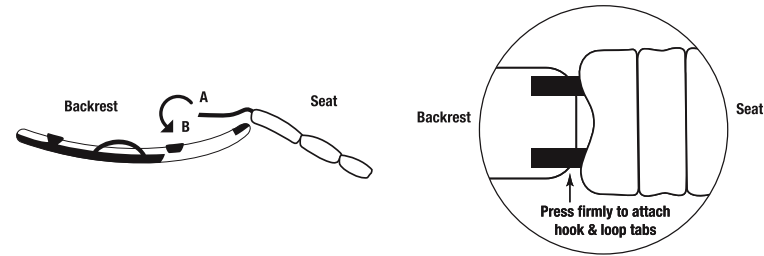

## Obus Ultra Forme™ Seat Use and Care Tips

- The one-size-fits-all Obus Ultra Forme™ Seat can fit any chair, anywhere.
- The seat should never overlap or fold over the front of any chair. If it does, simply remove the front zippered section.
- Sponge-wash the cover of your seat; DO NOT rub excessively.
- DO NOT place the seat or its cover in the washing machine.

## Materials

- 100% polycarbonate frame
- 100% polyurethane foam fill (CA 117 fire retardant)
- 100% polyester cover

## LIMITED ONE YEAR WARRANTY

Conair will repair or replace (at our option) your unit free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed below, together with your purchase receipt and \$10.00 for postage and handling. California residents need only provide proof of purchase and should call 1-800-366-0937 for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

**ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO, THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12-MONTH DURATION OF THIS WRITTEN, LIMITED WARRANTY.** Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

**IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights, which vary from state to state.

For information on any Conair product call: 1-800-3-CONAIR  
or visit us on the web at: [www.conair.com](http://www.conair.com)

©2007 Conair Corporation

**Service Center**  
Conair Corporation  
Service Department  
7475 N. Glen Harbor Blvd.  
Glendale, AZ 85307

07NP124

# OBUS<sup>ultra</sup>FORME™

SEAT

*Designed for back pain sufferers, including those with chronic back pain*



**OBUSFORME®**  
THE BODYCARE EXPERTS™  
BY **CONAIR™**

SE1

IB-7547

**INSTRUCTION BOOKLET FOR MODEL SE1**

For your safety and continued enjoyment of this product, always read the instruction book carefully before using.